

Evaluating various disintegrants regarding their performance in orally disintegrating tablet formulations

Thorsten Agnese¹, Florian Bang¹, Thorsten Cech¹

¹European Application Lab, Pharma Solutions, BASF SE, Ludwigshafen, Germany

Corresponding author: thorsten.cech@basf.com

INTRODUCTION

In regard to line extensions or improving administration convenience, orally disintegrating tablets (ODTs) have become a popular dosage form over the last years [1]. Even though formulators have some ready-to-use aids on hand, allowing quick and simple drug formulation [2], some active pharmaceutical ingredients raise the need to employ customised formulation development.

An ODT administration calls for fast disintegrating dosage forms offering a pleasant mouth sensation to gain patient's compliance. Hence, the selection of the disintegrant is crucial for the success of the final product. The aim of this study was to compare the performance of various disintegrants in a set ODT formulation [3].

MATERIALS AND METHODS

Lactose (GranuLac® 230, Meggle Pharma) was agglomerated (wet granulation) with native maize starch (C*PharmGel™, Cargill) to be used as filling material. The wet granulation process was conducted in a high shear mixer (Diosna P 1/6) applying an impeller speed of 200 rpm and a chopper speed of 2,000 rpm. The binder (2.0% w/w final granules) was added as aqueous paste within 120 s, followed by a granulation time of 180 s [3].

The wetted agglomerates were passed through an oscillating sieving machine (w=1.6 mm, AR400, ERWEKA), dried on a tray (ambient conditions), and finally passed through a sieve (w=0.8 mm) [3].

The following disintegrants were investigated with a concentration of 5.0% in the tableting blend: six grades of cross-linked poly(vinyl pyrrolidone) differing in particle size (Kollidon® CL, Kollidon® CL-F, Kollidon® CL-SF, Kollidon® CL-M, all BASF; Polyplasdone® XL, Polyplasdone® XL-10, both Ashland), croscarmellose sodium (Ac-Di-Sol®, FMC) and sodium starch glycolate (Explotab®, JRS).

Eventually, 0.5% magnesium stearate (Bärlocher) was added as lubricant to the tableting blend.

The compression was done using a single punch press XP 1 (Korsch) equipped with flat faced, faceted punches with a diameter of 8.0 mm. Compression forces of 2 to 8 kN were applied at a tableting speed of 20 tablets per minute.

The tablets were characterised (n=20) using a multi-tester (HT100, Sotax). The disintegration time (n=6) was determined (ERWEKA ZT 74) in demineralised water (37°C ±1 K).

RESULTS AND DISCUSSION

Generally, the characteristic of a tablet is strongly impacted by the disintegrant. Tensile strength and friability for instance (as essential features) are markedly influenced by particle size of the disintegrant chosen whereas smaller particles tend to lead to tablets of higher strength [4]. As regards to the ODT application, the texture of the disintegrated tablet presented to the tongue is worth noting [2].

In this regard, the latter characteristic is influenced by the wetted and swollen material whereas tensile strength is impacted by the particle size of the dry powder.

According to this, Kollidon® CL-M appears to be a promising candidate. But this grade is micronised and therefore lost most of its power to act as a strong disintegrant. Therefore, Kollidon® CL-SF was selected for this investigation (Figure 1).

As second cross-linked poly(vinyl pyrrolidone) grade Polyplasdone® XL-10 was chosen, since it also presented smaller particles compared to the regular Polyplasdone® XL grade. Furthermore, Ac-Di-Sol® and Explotab® were used, even though both products presented very uncomfortable sand in the mouth sensation due to the size of the swollen particles (Figure 2).

The visual appearance of all investigated disintegrants is shown in the SEM images Figure 3-10.

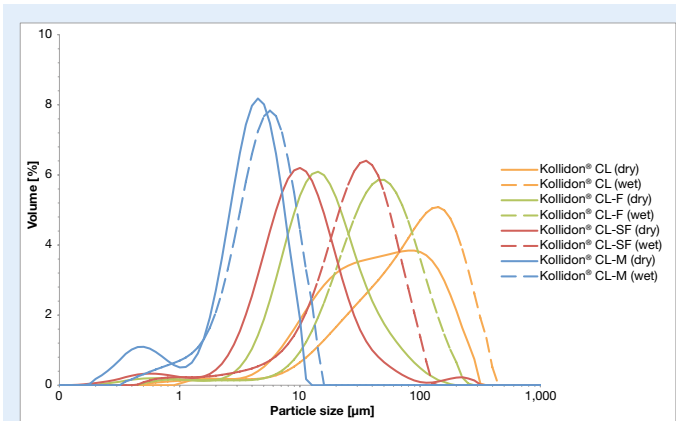


Figure 1.
Particle size distribution of the four different Kollidon® CL grades determined via laser diffraction in a dry and in a wet state (Mastersizer 2000, Malvern Instruments).

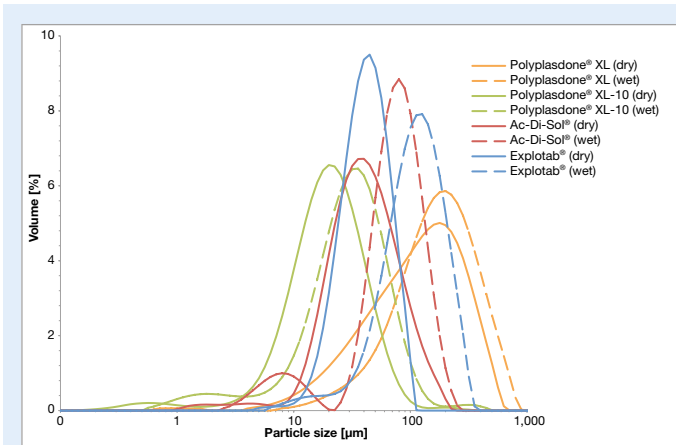


Figure 2.
Particle size distribution of four different disintegrants determined via laser diffraction in a dry and in a wet state (Mastersizer 2000, Malvern Instruments).

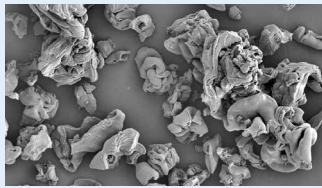


Figure 3.
Scanning electron microscopy (SEM) image of Kollidon® CL (SE, 5 kV).

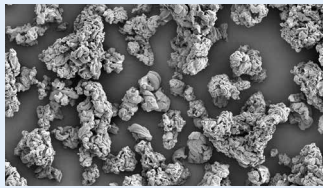


Figure 4.
Scanning electron microscopy (SEM) image of Kollidon® CL-F (SE, 5 kV).

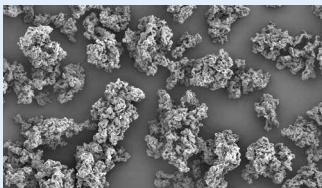


Figure 5.
Scanning electron microscopy (SEM) image of Kollidon® CL-SF (SE, 5 kV).

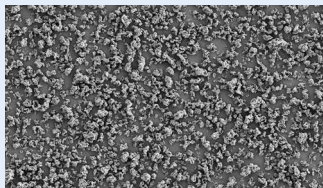


Figure 6.
Scanning electron microscopy (SEM) image of Kollidon® CL-M (SE, 5 kV).

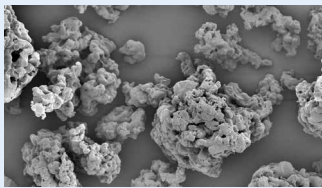


Figure 7.
Scanning electron microscopy (SEM) image of Polyplasdone® XL (SE, 5 kV).

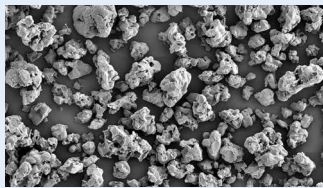


Figure 8.
Scanning electron microscopy (SEM) image of Polyplasdone® XL-10 (SE, 5 kV).

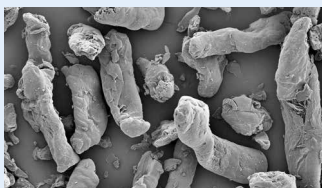


Figure 9.
Scanning electron microscopy (SEM) image of Ac-Di-Sol® (SE, 5 kV).

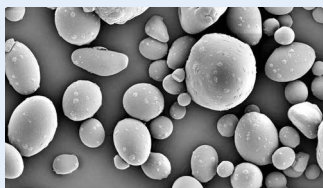


Figure 10.
Scanning electron microscopy (SEM) image of Explotab® (SE, 5 kV).

When compressing the different tableting blends to ODTs, a clear functionality of tensile strength on compression force could be found. Interestingly, adding a disintegrant to the tableting formulation hardly changed or even decreased the strength of the ODT. However, this was decisively different for Kollidon® CL-SF. Resulting tensile strength for formulations containing this particular excipient were distinctly higher compared to the disintegrant-free formulations; particularly at high compression forces (Figure 11).

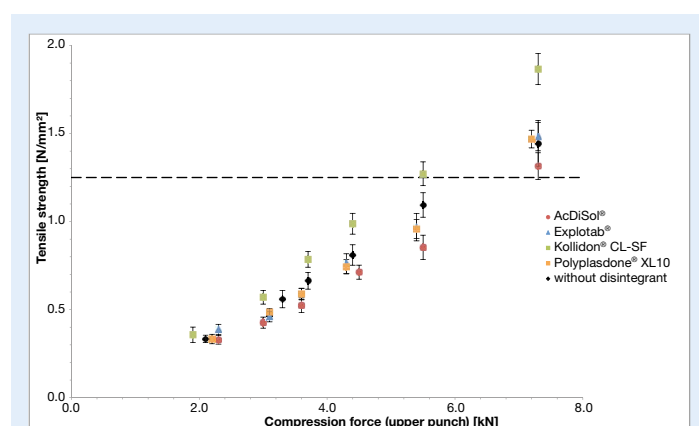


Figure 11.
Tensile strength of ODTs containing different disintegrants (or no disintegrant) as function of compression force (mean values (n=20), \pm SD).

Although tensile strength was found to be high, these tablets offered the shortest disintegration time of all formulations. And even more interesting, latter values were always found at merely about 10 s independent of the compression force applied (Figure 12).

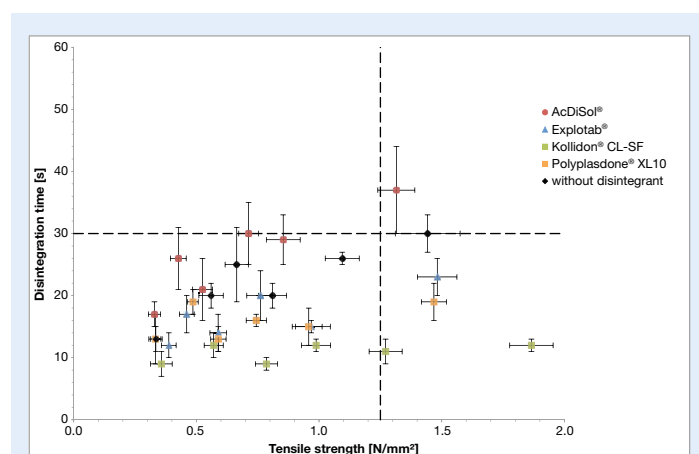


Figure 12.
Disintegration time of tablets, containing 10% of disintegrant (mean value (n=6), \pm SD).

CONCLUSION

Several advantages can be utilised when employing Kollidon® CL-SF as disintegrant in ODT tableting formulations: firstly, it presents a superior mouth sensation due to its small particle size, secondly, it leads to tablets of high tensile strength, and thirdly, it allows very fast disintegration independent of compression force applied.

REFERENCES

- [1] Bohnacker, R.; Streil, F.; Schweizer, S.; Müller, I.; Determination of the disintegration time of mouth melt tablets with texture analyser method; Pharm. Ind. 67 (3), 327-35 (2005).
- [2] Kruse, S.; Gebert, S.; Kolter, K.; et. al; Development of orally disintegrating tablets based on a new excipient, 2007 AAPS Annual Meeting and Exposition; November 11-15, 2007; San Diego (CA), USA
- [3] Agnese, Th., Cech, Th.; Evaluating various wet binders to gain lactose based agglomerates applicable for orally disintegrating tablet formulations; 1st European Conference on Pharmaceutics, April 13-14, 2015; Reims, France
- [4] Bang, F., Cech, Th.; Investigating the impact of particle size of different crospovidone grades on tablet characteristics; 9th PBP World Meeting, March 31 – April 3, 2014; Lisbon, Portugal

Kollidon® CL-SF

The finest grade of water-insoluble crospovidone disintegrant

- ✓ Particularly suitable for small and orally disintegrating tablets
- ✓ Creation of smooth tablet surfaces and cream-like mouth feel
- ✓ Suitable for wet granulation due to its strong ability to absorb water
- ✓ Compliant with current monographs for crospovidone (Type B in Ph. Eur. due to particle size)

Our service offer

We are providing in-depth expertise in all steps of the production of solid and liquid oral dosage forms.

The combination of our broad portfolio of functional excipients and our expert know-how enables you to create value-adding and unique formulations.

For more information visit us on
www.pharma.basf.com

For sample requests contact us at
pharma.solutions@basf.com

